



User Manual

QuantuMoonLight



Table of contents

1. Introduction	4
2. System description	4
3. System requirements	4
4. Overview	4
5. Registration	5
6. Authentication	5
7. User Profile Management	6
8. System functionalities	6
9. Dataset requirements	11
10. Community	11



List of Figures

Figure 1: The registration form of the platform.....	5
Figure 2: The authentication form of the platform.....	5
Figure 3: The homepage on QuantuMoonLight. The red arrow shows where the user area is.....	6
Figure 4: User profile in which the user can change the personal information, i.e., the IBM Token.....	6
Figure 5: Setup of the dataset for Quantum Computing.	7
Figure 6: List of classification and regression algorithms available in QuantuMoonLight.	8
Figure 7: Advanced setting options.....	8
Figure 8: The preprocessing steps.....	9
Figure 9: Feature extraction, Feature Selection and Prototype Selection.	10
Figure 10: Example of PegasosQSVC, Quantum Neural Network, SVC, and SVR Hyperparameters.....	11
Figure 11: The QuantuMoonLight community page.	12

1. Introduction

This manual represents a guide for the use of the system, making the comprehension of QuantuMoonLight easier and providing a preview of the functionalities.

2. System description

Quantum Machine Learning was born as a fusion of two ideas: Machine Learning and Quantum Computing. Despite the growing interest in the topic, few tools are useful to simplify the pipeline of quantum machine learning. Moreover, these solutions cannot be used by people without experience in software development. So, our work aims to offer a simple and intuitive product that can be used without deep programming knowledge.

3. System requirements

QuantuMoonLight is supported by any operating system: Windows, Linux, or Mac. The only requirement is a web browser in which the user can use QuantuMoonLight's website.

4. Overview

The platform required a free registration to use all the functionalities. You can find a short description of the system on the website's homepage.

The navbar on the top allows you to reach the services offered by our system:

- **Home:** contains a short description of the platform.
- **Quantum ML:** contains the main functionalities of the system.
- **Get Started:** contains information useful on the first tool usage.
- **Community:** allows you to interact with other users.
- **User Area:** contains your personal information and additional service, i.e., **Statistics and history of operations**, allowing you to see the history and the statistics on your experiments and compare them.

The footer contains information on the university in which the platform was developed, i.e., the University of Salerno. Finally, the page About Us contains information about the development team and the newsletter subscription form.

5. Registration

Regarding the registration on the platform, the user must select the User Area, click on “Register Here” and fill in all the requested fields. Please, consider inserting a valid e-mail and IBM token to use the system correctly. Please note that the password will be encrypted before being saved in the database (Figure 1).

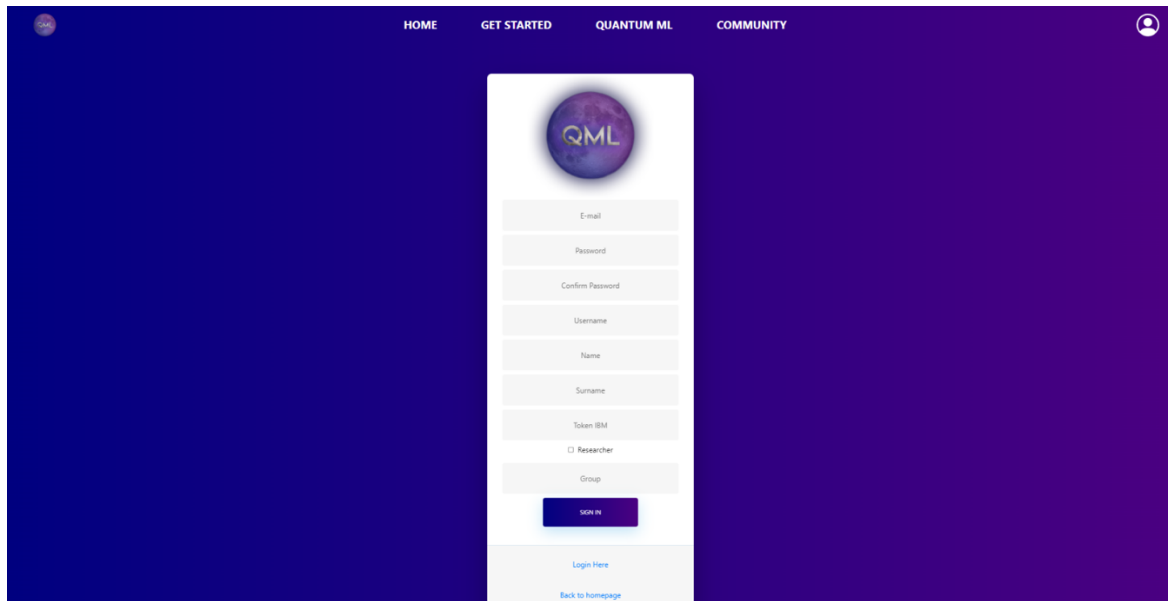
The image shows a registration form on a dark blue background. At the top, there is a navigation bar with links: HOME, GET STARTED, QUANTUM ML, and COMMUNITY. A user profile icon is in the top right corner. The form itself is a white box in the center, featuring the QML logo at the top. Below the logo are input fields for E-mail, Password, Confirm Password, Username, Name, Surname, and Token IBM. There is a checkbox for 'Researcher' and a 'Group' field. A blue 'SIGN IN' button is positioned below these fields. At the bottom of the form, there are links for 'Login Here' and 'Back to homepage'.

Figure 1: The registration form of the platform.

6. Authentication

The access to the platform is through the User Area by entering username and password (Figure 2).

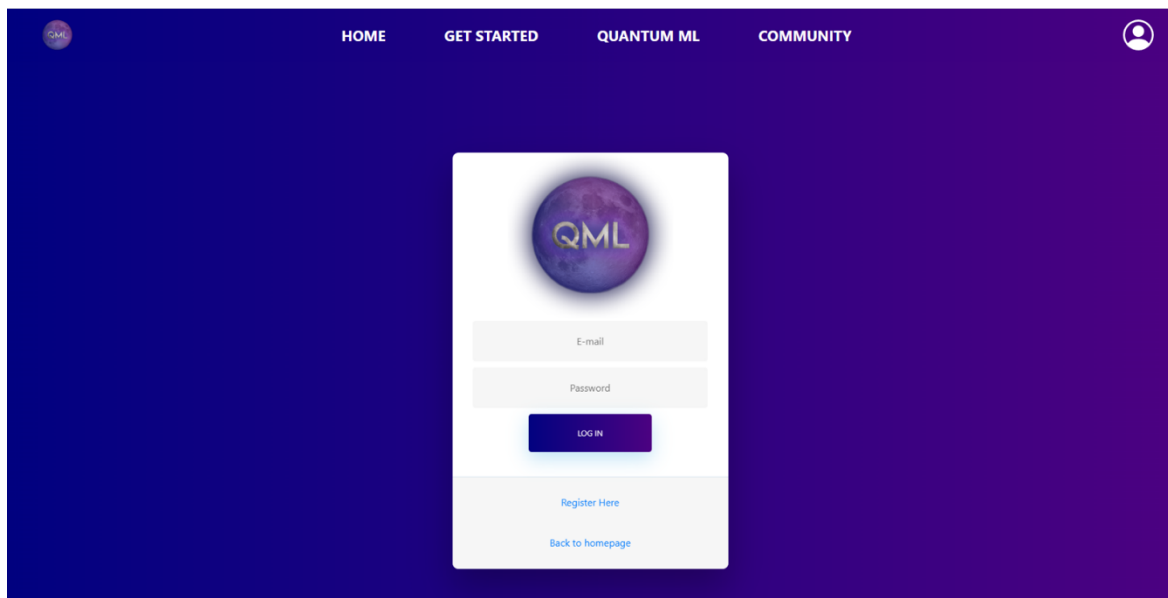
The image shows an authentication form on a dark blue background. At the top, there is a navigation bar with links: HOME, GET STARTED, QUANTUM ML, and COMMUNITY. A user profile icon is in the top right corner. The form is a white box in the center, featuring the QML logo at the top. Below the logo are input fields for E-mail and Password. A blue 'LOG IN' button is positioned below these fields. At the bottom of the form, there are links for 'Register Here' and 'Back to homepage'.

Figure 2: The authentication form of the platform.

To logout the user must go to the User Area and click the appropriate “Log-Out” button.

7. User Profile Management

To change your profile information, you need to go in the user area, as showed in Figure 3. Once the user goes in the user area, Figure 4 shows the user profile in which the user can see and change the personal information.

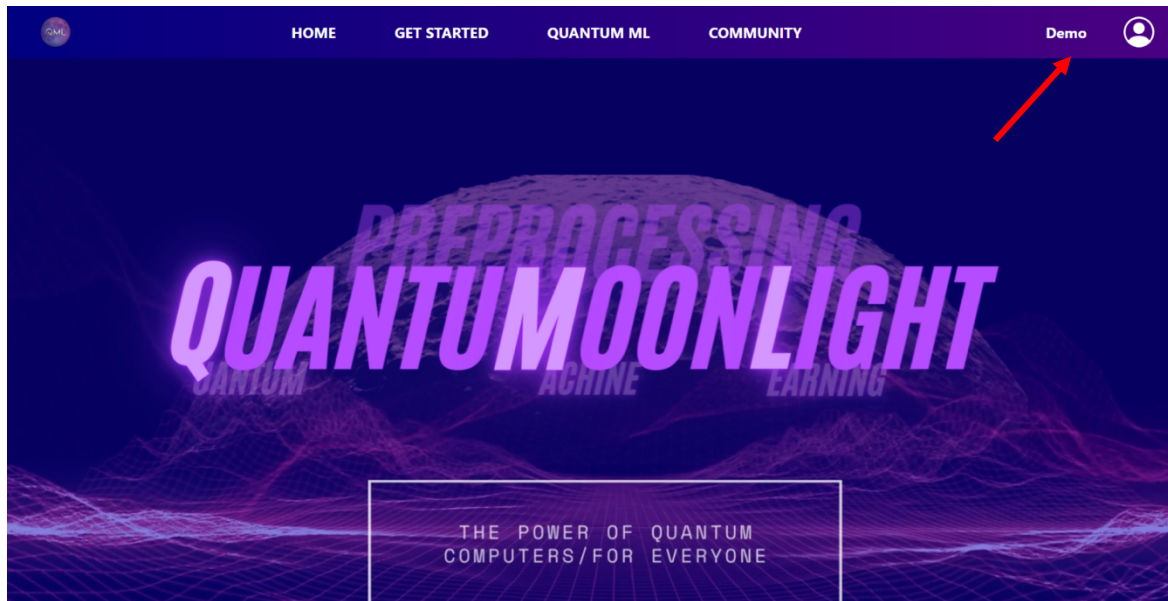


Figure 3: The homepage on QuantuMoonLight. The red arrow shows where the user area is.

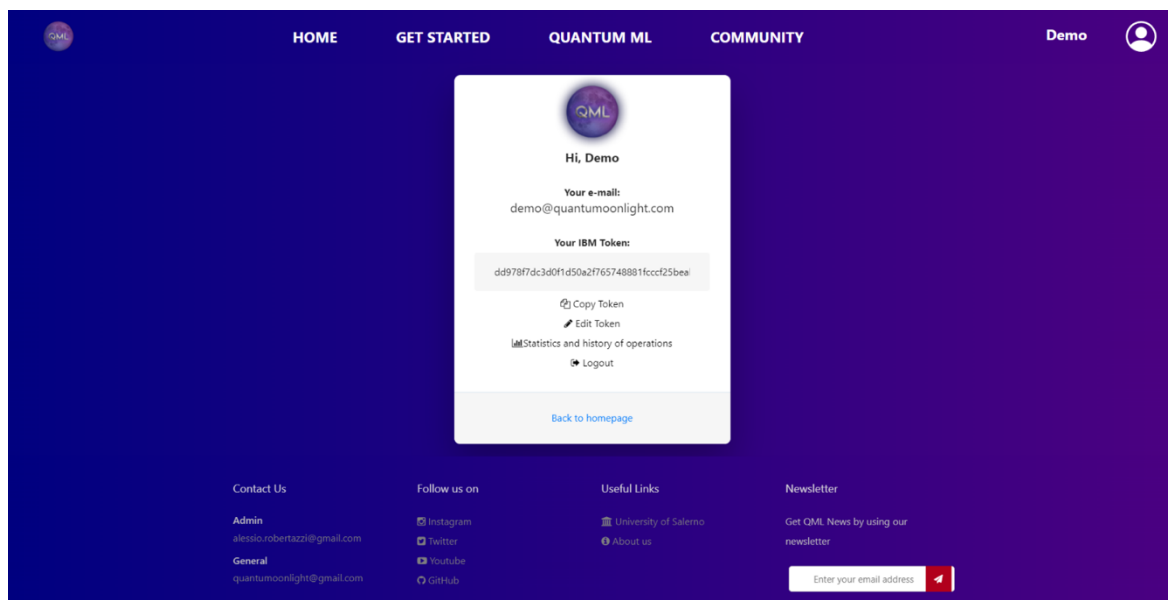


Figure 4: User profile in which the user can change the personal information, i.e., the IBM Token.

8. System functionalities

The Quantum ML page contains information through which the user can decide on validation, preprocessing and classification steps. As the reader may notice in Figure 5, there is a box where the

user can load Training (with the possibility of also selecting an automatic split), Test, and Prediction sets. It is also possible to preprocess data to reduce columns through the Feature Extraction by PCA or the Feature Selection by K-best. Finally, it is also possible to minimize rows with the Prototype Selection.

Load your Dataset for Quantum Classification

Training Set to upload [Choose File](#)

Test Set to Upload [Choose File](#)

Prediction Set to Upload [Choose File](#)

☐ Automatically Split the Training Set ⓘ

☐ Reduce columns with Feature Extraction(PCA)

☐ Reduce columns with Feature Selection(K Best) ⓘ

☐ Reduce rows with Prototype Selection ⓘ

Figure 5: Setup of the dataset for Quantum Computing.

Figure 6 shows a list of different algorithms implemented in our tool. In particular, there are QSVM, QSVC, PegasosQSVC, and QNN for quantum classification. Then, QSVR, and VQR for quantum regression; As the reader may see SVC, KNN, Decision Tree, Naive Bayes, and Random Forest for canonical classification. Finally, SVR and Linear Regression for canonical regression.

Load your Dataset for Quantum Classification

Training Set to upload

Test Set to Upload

Prediction Set to Upload

☐ Automatic

☐ Reduce columns

☐ Reduce columns

☐ Reduce rows

Select the algorithm:

- None
- Quantum Classifiers**
- QSVM
- QSVC
- Pegasos QSVC
- Quantum Neural Network
- Quantum Regressors**
- QSVR
- VQR
- Canonical Classifiers**
- SVC
- K Neighbors
- Decision Tree
- Naive Bayes
- Random Forest
- Canonical Regressors**
- SVR
- Linear Regression

Figure 6: List of classification and regression algorithms available in QuantuMoonLight.

Users can set advanced options such as K-Fold Cross Validation by choosing the “K” parameter (default is 10) and entering the number of columns and rows to be reduced, as shown in Figure 7.

Validation

☐ Simple Train-Test Split

☒ K-fold Cross Validation

Insert k parameter:

Figure 7: Advanced setting options.

Figure 8 shows the preprocessing process; the reader may notice that the users can select Data Balancing with SMOTE and/or Data Imputation and choose between Min-Max Scaling and Standard Scaling.

The image shows a software interface for 'Advanced Options'. At the top is a blue button with a wrench icon and the text 'Advanced Options'. Below this is a large, light-gray rectangular area containing two sections. The first section is titled 'Validation' in purple. It contains two radio button options: 'Simple Train-Test Split' (which is selected) and 'K-fold Cross Validation'. Each option has a small circular information icon (i) to its right. The second section is titled 'Preprocessing' in purple. It contains four checkbox and radio button options: 'Data Imputation' (checkbox), 'Data Balancing' (checkbox), 'MinMax Scaling' (radio), and 'Standard Scaling' (radio). Each option also has a small circular information icon (i) to its right.

Figure 8: The preprocessing steps.

Moreover, it is possible to choose also the number of columns for Feature Extraction and Feature Selection and the number of rows for Prototype Selection to reduce (Figure 9).

The screenshot shows a vertical stack of three configuration panels. The top panel, titled 'Feature Extraction', contains the text 'Insert nr of columns to reduce Feature Extraction:' followed by a text input field containing the number '2' and an information icon. The middle panel, titled 'Feature Selection', contains the text 'Insert nr of columns to reduce Feature Selection:' followed by a text input field containing the number '2' and an information icon. The bottom panel, titled 'Prototype Selection', contains the text 'Insert nr of rows to reduce:' followed by a text input field containing the number '100' and an information icon.

Figura 9: Feature extraction, Feature Selection and Prototype Selection.

Finally, for Classification and Regression algorithms, the users may change the hyperparameters and choose the Backend of IBMQ family (Figure 10). For this option, you will need to see the backend associated with the token. These backends are divided into Simulator and Quantum System. After uploading the dataset, the system will redirect the user to a page where they can download validation and preprocessing results. If you have selected an algorithm for classification/regression, you will receive the results via mail.

The screenshot displays two side-by-side configuration panels. The left panel, titled 'Hyperparameter', is for 'Pegasos QSVN' and includes two settings: 'Insert PegasosQSVN C:' with a text input field containing '1000' and an information icon, and 'Insert PegasosQSVN tau:' with a text input field containing '100' and an information icon. The right panel, also titled 'Hyperparameter', is for 'Quantum Neural Network' and includes three settings: 'Insert NN optimizer:' with a dropdown menu showing 'SLSQP' and an information icon, 'Insert NN loss function:' with a dropdown menu showing 'Squared Error' and an information icon, and 'Insert max iterations optimizer:' with a text input field containing '100' and an information icon.

Hyperparameter

Insert SVC kernel: ⓘ

Insert SVC C: ⓘ

Hyperparameter

Insert SVR kernel: ⓘ

Insert SVR C: ⓘ

Figure 10: Example of PegasosQSVC, Quantum Neural Network, SVC, and SVR Hyperparameters.

9. Dataset requirements

The system does not check the quality of your dataset and does not provide any data cleaning operations. Therefore, the user who will upload the datasets, should take a look at the following requirements:

- The dataset should be a CSV file.
- The dataset should contain only numbers and, at most, the “NaN” value, separated by commas.
- The dataset does not include the column id and the name of the features.
- Each row must have the same number of values.
- The dataset should not have strings or empty strings.
- The dataset useful to predict must be coherent at the training set, i.e., the same number of the features.
- The dataset useful to predict must contain just one column less than the training and test datasets.

10. Community

In the last section of our system, i.e., the community, the registered users could interact with each other by joining a discussion and/or sharing their results. To create a new thread, you must login to

the system. Then, as shown in Figure 11, the user could choose to share an article or an experiment and /or comment the other posts.

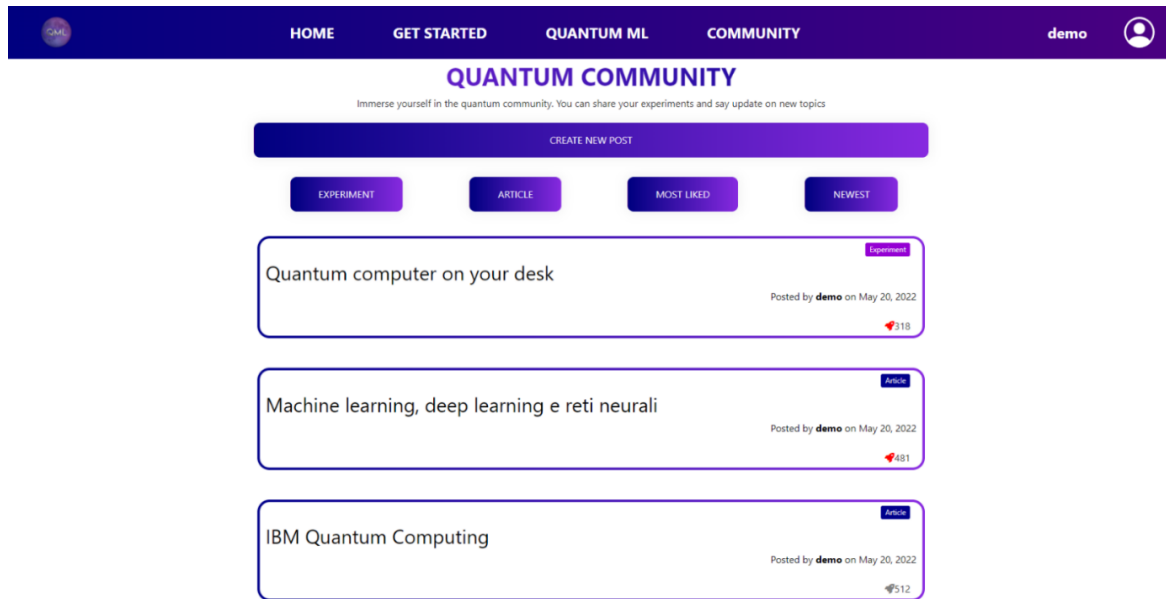


Figure 11: The QuantuMoonLight community page.